

Slant Back Whirlpools

Division of Acorn Engineering Company











SB Series 100 Gallon Slant Back Mobile Whirlpool

SB-100-M



SB-100-M

Whitehall SB-100-M Slant Back Mobile Whirlpool is recommended for treatment of extremities, hip and back area. A slant back design allows for the user to comfortably recline while soaking in the tank.

Tank Capacity is 100 gallons.

Whirlpool Tank is 76"L x 40"W x 24"D and is seamless welded construction, fabricated from heavy gauge, type 304 stainless steel, polished to a satin finish. The tank features a slant back to allow the user to recline while soaking in the tank. A 3" wide rim allows the user to also rest their arms while soaking. Tank bottom is seamless, coved (rounded) design for effective cleaing to minimize bacteria build-up.

Turbine Raising and Lowering Device permits adjustment of desired height and direction of water agitation towards the area of the body to be treated. It is spring balanced and provided with a locking device.

Turbine Assembly functions both as agitator and emptying device. It is equipped with Underwriters Laboratory listed 1/2 HP pump motor protected with an ABS plastic shell, automatic thermal overload protector, and lifetime-sealed bearings. All parts contacting the water are stainless steel, chrome plated brass or bronze.

Drain option -SDP has a bottom drain, leading to a separate, concealed water evacuation pump. A hose is attached to the side of the unit to facilitate the evacuation.

Thermometer is provided with a dual scale, adjustable calibration. Stem inserts into a retainer to prevent rattle during whirlpool agitation.

GUIDE SPECIFICATIONS

Provide Whitehall Mobile Whirlpool (specify model number and options) with a 100-gallon capacity, 76"L x 40"W x 24"D. Whirlpool shall be fabricated from heavy gauge, type 304 stainless steel. Construction shall be seamless welded and exposed surfaces shall be polished to satin finish. The tank features a slant back to allow the user to recline while soaking in the tank. The 3" rim allows the user to rest their arms while soaking. Bottom of tank shall be a seamless, rounded coved design to minimize bacteria build-up. Provide a turbine assembly with raising and lowering device which functions both as agitator and emptying device, permits adjustment of desired height, direction of water agitation, and can be locked into place. Turbine assembly shall be UL listed 1/2 HP pump motor with automatic thermal overload protector and lifetime-sealed bearings. Provide a 3" dial thermometer with a stem retainer to prevent rattle when whirlpool is in operation.

Please visit www.whitehallmfg.com for most current specifications.





Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 782-7706 • (626) 968-6681 • Fax (626) 855-4862

Updated: 10/08/18



Slant Back Whirlpools

MODEL NUMBER AND OPTION SELECTIONS

Series

✓ SB Slant Back Series

Capacity

₫ -100 Gallons

✓ -M Mobile Whirlpool

Turbine Electrical Requirements

GFI Receptacle

☐ 115VAC / 50Hz ☐ 115VAC / 60HZ □ 230VAC / 50Hz ☐ 230VAC / 60HZ Country of Use

For Mixing Valve:

Volts: 115

Refer to supplemental MXT15 sheet for details.

Suffix Options (See Whirlpool Options Section for Details)

Auxiliary Drain Hose Assembly □ -ADH

☐ -ALC Adjustable Low Chair □ -ALS Adjustable Low Stool

□ -EG Enviro Glaze Color Powder Coating -

Specify Color:_

(Special Finishes Consult Factory)


☐ -GFI Ground Fault Interrupter ☐ -SDP Separate Drain Pump

ELECTRIC REQUIREMENTS: Amps: 6.9/6.2 Hertz: 50/60

ELECTRIC REQUIREMENTS: Amps: 3.5/3.1 Hertz: 50/60 Volts: 230 GFI Receptacle Certified per US & Canadian Medical Standards: UL 60601-1 & CSA C22.2 No


European Standards; (E IEC 60601-1




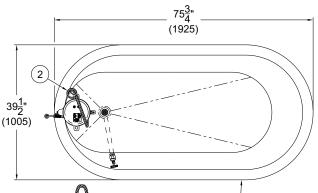

NOTES:

- 1. Slant Back Whirlpool
- 2. Thermometer
- 3. 1/2" HP Pump Motor
- 4. Turbine Raising & Lowering Device
- 5. Auxiliary Drain

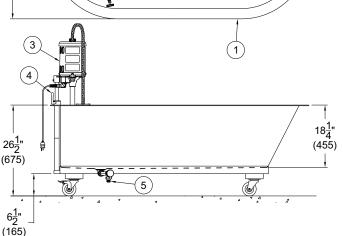

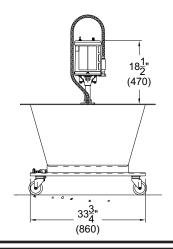

▲ WARNING: Cancer and Reproductive Harm - www.P65Warnings.ca.gov

All dimensions are nominal and subject to manufacturer's change without notice. Whitehall assumes no responsibility fo use of void or superseded data © Whitehall Manufacturing, Member of Morris Group International.

Please visit www.whitehallmfg.com for most current specifications

SELECTION SUMMARY & APPROVAL FOR MANUFACTURING			
Model Number & Options		Quantity	
Company		Date	
Contact	Title		
Approval for Manufacturing/Signature			

Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 782-7706 • (626) 968-6681 • Fax (626) 855-4862

Updated: 10/08/18