

Lo-Boy Whirlpools

Division of Acorn Engineering Company



A MEMBER OF





L Series 75 Gallon Mobile Whirlpool

L-75-M



L-75-M

Whitehall L-75-M Mobile Whirlpools are recommended for treatment of extremities, hip and back area. It is ideal for pediatric and geriatric applications imparting a greater degree of patient comfort and safety.

Tank Capacity is 75 gallons.

Whirlpool Tank is 52"L x 24"W x 18"D and is seam-less welded construction, fabricated from heavy gauge, type 304 stainless steel, polished to a satin finish. Tank rim is reinforced with a concealed stainless steel rod. Tank bottom is a seamless, coved (rounded) design for effective cleaning to minimize bacteria buildup.

Heavy Duty Swivel Casters support all mobile models. Rear pair of casters are lockable for stabilizing whirlpool during treatment.

Turbine Raising and Lowering Device permits adjustment of desired height and direction of water agitation towards the area of the body to be treated. It is spring balanced and provided with a locking device.

Turbine Assembly functions both as agitator and emptying device. It is equipped with Underwriters Laboratory listed 1/2 HP jet pump motor protected with an ABS plastic shell, automatic thermal overload protector, and lifetime-sealed bearings. All parts contacting the water are stainless steel, chrome plated brass or bronze.

Drain option -SDP has a bottom drain, leading to a separate, concealed water evacuation pump. A hose is attached to the side of the unit to facilitate the evacuation.

Thermometer is provided with a dual scale, adjustable calibration, and a stem inserts into a retainer to prevent rattle during whirlpool agitation.

GUIDE SPECIFICATIONS

Provide Whitehall Mobile Whirlpool (specify model number and options) with a 75-gallon capacity, 52"L x 24"W x 88"D. Whirlpool shall be fabricated from heavy gauge, type 304 stainless steel. Construction shall be seamless welded and exposed surfaces shall be polished to satin finish. Bottom of tank shall be a seamless, rounded coved design to minimize bacteria build up. Whirlpool shall have heavy duty swivel casters with locking device on rear pair of casters. Provide a turbine assembly with raising and lowering device which functions both as agitator and emptying device, permits adjustment of desired height, direction of water agitation, and can be locked into place. Turbine assembly shall be UL listed 1/2 HP jet pump motor with automatic thermal overload protector and lifetime-sealed bearings. Provide a thermometer with a stem retainer to prevent rattling when whirlpool is in operation.

Please visit www.whitehallmfg.com for most current specifications.





Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 782-7706 • (626) 968-6681 • Fax (626) 855-4862

Updated: 10/03/18



Lo-Boy Whirlpools

MODEL NUMBER AND OPTION SELECTIONS Series ✓ L Lo-Boy Series **Capacity** ₫ -75 Gallons Type ✓ -M Mobile Whirlpool **Turbine Electrical Requirements** ☐ 115VAC / 50Hz ☐ 115VAC / 60HZ ☐ 230VAC / 50Hz ☐ 230VAC / 60HZ


Suffix O	ptions (See Whirlpool Options Section for Details)		
🔲 -ADH	Auxiliary Drain Hose Assembly		
-ALC	Adjustable Low Chair		
-ALS	Adjustable Low Stool		
🔲 -ARS	Arm Rest Support		
☐ -AS3	Adjustable Suspension Seat		
-CTS3	Combination Table with Side Seats		
🔲 -EG	Enviro Glaze Color Powder Coating -		
	Specify Color:		
	(Special Finishes Consult Factory)		
🔲 -GFI	Ground Fault Interrupter		
🔲 -HR	Head Rest, Adjustable		
-PCS	Pump Strainer Cover		
-PHS2	Patient Handling System (Chair Type)		
🔲 -RHG	Rubber Hand Grips		
🔲 -SDP	Separate Drain Pump		

☐ -TTS2 Tank Top Seat

For Mixing Valve:

Country of Use

Refer to supplemental MXT15 sheet for details.

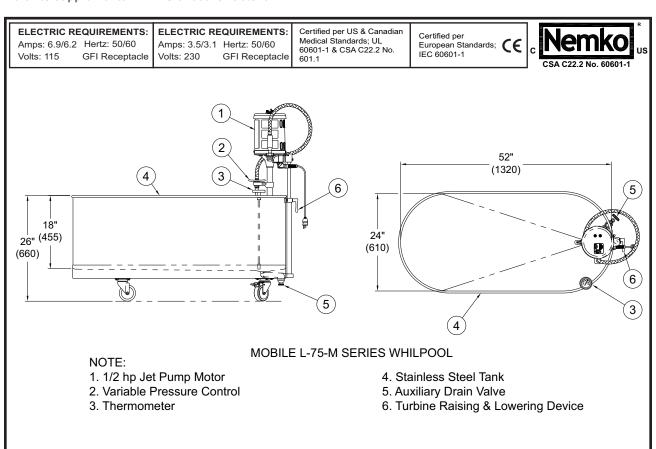

All dimensions are nominal and subject to manufacturer's change without notice. Whitehall assumes no responsibility fo use of void or superseded data © Whitehall Manufacturing, Member of Morris Group International.

Please visit www.whitehallmfg.com for most current specifications

SELECTION SUMMARY & APPROVAL FOR MANUFACTURING			
Model Number & Options		Quantity	
Company		Date	
Contact	Title		
Approval for Manufacturing/Signature			

Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 782-7706 • (626) 968-6681 • Fax (626) 855-4862

Updated: 10/03/18