

Hi-Boy Whirlpools

MORRIS GROUP

Use



H Series 60 Gallon Mobile Whirlpool

H-60-M



Whitehall H-60-M Mobile Whirlpools are recommended for shallow or deep body immersion offering complete coverage of the lower extremities as well as the lower back and hip.

Tank Capacity is 60 gallons.

Whirlpool Tank is 36"L x 20"W x 28"D and is seam-less welded construction, fabricated from heavy gauge, type 304 stainless steel, polished to a satin finish. Tank rim is reinforced with a concealed stainless steel rod. Tank bottom is a seamless, coved (rounded) design for effective cleaning to minimize bacteria buildup.

Heavy Duty Swivel Casters support all mobile models. Rear pair of casters are lockable for stabilizing whirlpool during treatment.

Turbine Raising and Lowering Device permits adjustment of desired height and direction of water agitation towards the area of the body to be treated. It is spring balanced and provided with a locking device.

Turbine Assembly functions both as agitator and emptying device. It is equipped with Underwriters Laboratory listed 1/2 HP jet pump motor protected with an ABS plastic shell, automatic thermal overload protector, and lifetime-sealed bearings. All parts contacting the water are stainless steel, chrome plated brass or bronze.

Drain option -SDP has a bottom drain, leading to a separate, concealed water evacuation pump. A hose is attached to the side of the unit to facilitate the evacuation.

Thermometer is provided with a dual scale, adjustable calibration, and a stem inserts into a retainer to prevent rattle during whirlpool agitation.

GUIDE SPECIFICATIONS

Provide Whitehall Mobile Whirlpool (specify model number and options) with a 60-gallon capacity, 36"L x 20"W x 28"D. Whirlpool shall be fabricated from heavy gauge, type 304 stainless steel. Construction shall be seamless welded and exposed surfaces shall be polished to satin finish. Bottom of tank shall be a seamless, rounded coved design to minimize bacteria build up. Whirlpool shall have heavy duty swivel casters with locking device on rear pair of casters. Provide a turbine assembly with raising and lowering device which functions both as agitator and emptying device, permits adjustment of desired height, direction of water agitation, and can be locked into place. Turbine assembly shall be UL listed 1/2 HP jet pump motor with automatic thermal overload protector and lifetime-sealed bearings. Provide a thermometer with a stem retainer to prevent rattling when whirlpool is in operation.

Please visit www.whitehallmfg.com for most current specifications.





Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 782-7706 • (626) 968-6681 • Fax (626) 855-4862

Updated: 10/03/18



Hi-Boy Whirlpools

MODEL NUMBER AND OPTION SELECTIONS Series ☑ H Hi-Boy Series Capacity ☑ -60 Gallons Type ☑ -M Mobile Whirlpool Turbine Electrical Requirements □ 115VAC / 50Hz □ 115VAC / 60HZ □ 230VAC / 50Hz □ 230VAC / 60HZ

Refer to supplemental MXT15 sheet for details

Suffix Options (See Whirlpool Options Section for D		
	🔲 -ADH	Auxiliary Drain Hose Assembly
	-AHC1	Adjustable High Chair with Casters
	-AHC2	Adjustable High Chair with Rubber Tips
	🔲 -ARS	Arm Rest Support
	🔲 -AS1	Adjustable Suspension Seat
	-CTS1	Combination Table with Side Seats
	🔲 -EG	Enviro Glaze Color Powder Coating -
		Specify Color:
		(Special Finishes Consult Factory)
	🔲 -GFI	Ground Fault Interrupter
	-PCS	Pump Strainer Cover
	-PHS2	Patient Handling System (Chair Type)
	🔲 -RHG	Rubber Hand Grips
	-SDP	Separate Drain Pump
	-TTS1	Tank Top Seat

ELECTRIC REQUIREMENTS:
Amps: 6.9/6.2 Hertz: 50/60
Volts: 115 GFI Receptacle

Country of Use


For Mixing Valve:


ELECTRIC REQUIREMENTS:
Amps: 3.5/3.1 Hertz: 50/60
Volts: 230 GFI Receptacle

Certified per US & Canadian Medical Standards; UL 60601-1 & CSA C22.2 No.

Certified per European Standards; **CE** IEC 60601-1



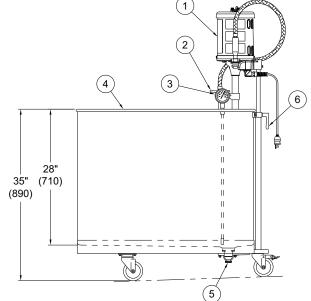

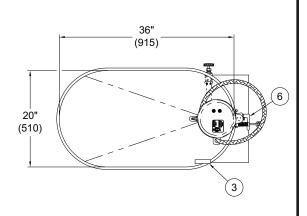

NOTES:

- 1. 1/2 hp Jet Pump Motor
- 2. Variable Pressure Control
- 3. Thermometer

- 4. Stainless Steel Tank
- 5. Auxiliary Drain Valve
- 6. Turbine Raising & Lowering Device

All dimensions are nominal and subject to manufacturer's change without notice. Whitehall assumes no responsibility fo use of void or superseded data © Whitehall Manufacturing, Member of Morris Group International.

Please visit **www.whitehallmfg.com** for most current specifications.

▲ WAI	RNING: Cancer and Reproductive Harm - www.P6	5Warnings.ca.gov			
SELECTION SUMMARY & APPROVAL FOR MANUFACTURING					
Model Number & OptionsQuantity					
Company	Date				
Contact	Title				
Approval for Manufacturing/Signature					

Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 782-7706 • (626) 968-6681 • Fax (626) 855-4862

MOBILE H-60-M SERIES WHIRLPOOL