

Hi-Boy Whirlpools

Division of Acorn Engineering Company











H Series 105 Gallon Mobile Whirlpool

H-105-M



H-105-M

Whitehall H-105-M Mobile Whirlpools are recommended for shallow or deep body immersion offering complete coverage of extremities as well as the entire back and shoulder area with greater patient comfort.

Tank Capacity is 105 gallons.

Whirlpool Tank is 48"L x 24"W x 28"D and is seam-less welded construction, fabricated from heavy gauge, type 304 stainless steel, polished to a satin finish. Tank rim is reinforced with a concealed stainless steel rod. Tank bottom is a seamless, coved (rounded) design for effective cleaning to minimize bacteria buildup.

Heavy Duty Swivel Casters support all mobile models. Rear pair of casters are lockable for stabilizing whirlpool during treatment.

Turbine Raising and Lowering Device permits adjustment of desired height and direction of water agitation towards the area of the body to be treated. It is spring balanced and provided with a locking device.

Turbine Assembly functions both as agitator and emptying device. It is equipped with Underwriters Laboratory listed 1/2 HP jet pump motor protected with an ABS plastic shell, automatic thermal overload protector, and lifetime-sealed bearings. All parts contacting the water are stainless steel, chrome plated brass or bronze.

Drain option -SDP has a bottom drain, leading to a separate, concealed water evacuation pump. A hose is attached to the side of the unit to facilitate the evacuation.

Thermometer is provided with a dual scale, adjustable calibration, and a stem inserts into a retainer to prevent rattle during whirlpool agitation.

GUIDE SPECIFICATIONS

Provide Whitehall Mobile Whirlpool (specify model number and options) with a 105-gallon capacity, 48"L x 24"W x 28"D. Whirlpool shall be fabricated from heavy gauge, type 304 stainless steel. Construction shall be seamless welded and exposed surfaces shall be polished to satin finish. Bottom of tank shall be a seamless, rounded coved design to minimize bacteria build up. Whirlpool shall have heavy duty swivel casters with locking device on rear pair of casters. Provide a turbine assembly with raising and lowering device which functions both as agitator and emptying device, permits adjustment of desired height, direction of water agitation, and can be locked into place. Turbine assembly shall be UL listed 1/2 HP jet pump motor with automatic thermal overload protector and lifetime-sealed bearings. Provide a thermometer with a stem retainer to prevent rattling when whirlpool is in operation.

Please visit **www.whitehallmfg.com** for most current specifications.






Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 782-7706 • (626) 968-6681 • Fax (626) 855-4862



Hi-Boy Whirlpools

MODEL NUMBER AND OPTION SELECTIONS Suffix Options (See Whirlpool Options Section for Details) □ -ADH Auxiliary Drain Hose Assembly **Series** -AHC1 Adjustable High Chair with Casters H Hi-Boy Series ☐ -AHC2 Adjustable High Chair with Rubber Tips □ -ARS Arm Rest Support **Capacity** Adjustable Suspension Seat ☐ -AS2 ₫ -105 Gallons -CTS2 Combination Table with Seats **Type** Enviro Glaze Color Powder Coating -☐ -EG ✓ -M Mobile Whirlpool Specify Color:_ (Special Finishes Consult Factory) **Turbine Electrical Requirements** Ground Fault Interrupter ☐ -GFI ☐ 115VAC / 50Hz ☐ 115VAC / 60HZ Pump Strainer Cover -PCS ☐ 230VAC / 50Hz ☐ 230VAC / 60HZ -PHS2 Patient Handling System (Chair Type) □ -RHG Rubber Hand Grips Country of Use ☐ -SDP Separate Drain Pump □ -TTS2 Tank Top Seat For Mixing Valve: Refer to supplemental MXT15 sheet for details.

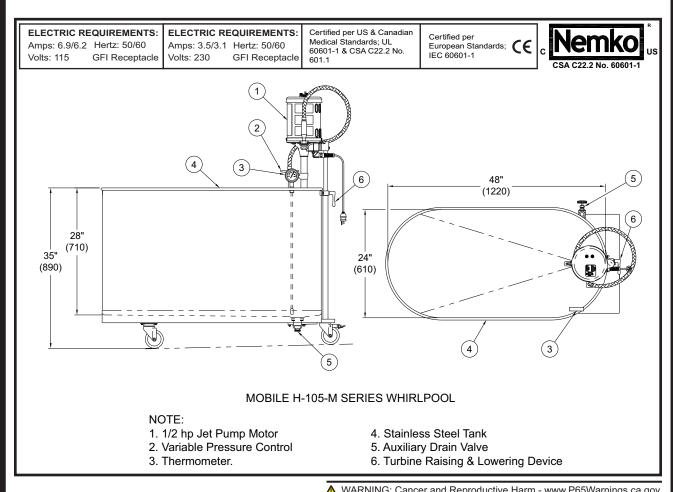

All dimensions are nominal and subject to manufacturer's change without notice. Whitehall assumes no responsibility fo use of void or superseded data © Whitehall Manufacturing, Member of Morris Group International.

Please visit **www.whitehallmfg.com** for most current specifications.

	A William C. Cancor and Represente Flan	www corranningc.ca.go.
SELECTION SUMMARY & APPROVAL FOR MANUFACTURING		
Model Nu	mber & Options	Quantity
Company		Date
Contact_	Title	
Approval for Manufacturing/Signature		

Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 782-7706 • (626) 968-6681 • Fax (626) 855-4862