

Extremity Whirlpools









E Series 22 Gallon **Stationary Whirlpool** with Pedestal

E-22-SP



E-22-SP

Whitehall E-22-S Stationary Whirlpool with **Pedestal** is recommended for treatment of the foot, ankle and leg, to below the knee, with better coverage and greater patient comfort.

Tank Capacity is 22 gallons.

Whirlpool Tank is 28"L x 15"W x 18"D and is seamless welded construction, fabricated from heavy gauge, type 304 stainless steel, polished to a satin finish. Tank rim is reinforced with a concealed stainless steel rod. Tank bottom is seamless, coved (rounded) design for effective cleaning to minimize bacteria buildup. Tank is fitted with a fixed raised pedestal for mounting

A 2" Combination Drain and Overflow, and a Filler Spout are provided.

Turbine Raising and Lowering Device permits adjustment of desired height and direction of water agitation towards the area of the body to be treated. It is spring balanced and provided with a locking device.

Turbine Assembly is equipped with Underwriters Laboratory listed 1/2 HP jet pump motor protected with an ABS plastic shell, and automatic thermal overload protector and lifetime-sealed bearings. All parts contacting the water are stainless steel, chrome plated brass or bronze.

Thermometer is provided with a dual scale, adjustable calibration, and a stem inserts into a retainer to prevent rattle during whirlpool agitation.

GUIDE SPECIFICATIONS

Provide Whitehall Stationary Whirlpool (specify model number and options) with a 22-gallon capacity, 28"L x 15"W x 18"D. Whirlpool shall be fabricated from heavy gauge, type 304 stainless steel. Construction shall be seamless welded and exposed surfaces shall be polished to satin finish. Tank shall have a reinforced rim and shall be fitted with a fixed raised pedestal. Bottom of tank shall be a seamless, rounded coved design to minimize bacteria build up. Provide a 2" combination drain and overflow, and a filler spout. Provide a turbine assembly with raising and lowering device which functions as an agitator and permits adjustment of desired height, direction of water agitation, and can be locked into place. Turbine assembly shall be UL listed 1/2 HP jet pump motor with automatic thermal overload protector and lifetime-sealed bearings. Provide a thermometer with a stem retainer to prevent rattle when whirlpool is in operation.

Please visit www.whitehallmfg.com for most current specifications.





Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 782-7706 • (626) 968-6681 • Fax (626) 855-4862

Updated: 10/02/18



Extremity Whirlpools

MODEL NUMBER AND OPTION SELECTIONS Series ☑ E Extremity Series **Capacity** ✓ -22 Gallons **Type**

✓ -SP Stationary Whirlpool with Pedestal **Turbine Electrical Requirements**

☐ 115VAC / 50Hz ☐ 115VAC / 60HZ ☐ 230VAC / 50Hz ☐ 230VAC / 60HZ

Country of Use


For Mixing Valve:

ELECTRIC REQUIREMENTS:

Refer to supplemental MXT15 sheet for details.

Suffix Options (See Whirlpool Options Section for Details) □ -AHC1 Adjustable High Chair with Casters ☐ -AHC2 Adjustable High Chair with Rubber Tips □ -ARS Arm Rest Support □ -EG Enviro Glaze Color Powder Coating -Specify Color: (Special Finishes Consult Factory) ☐ -GFI Ground Fault Interrupter □ -PCS Pump Strainer Cover

☐-SDO Separate Drain and Overflow Assembly

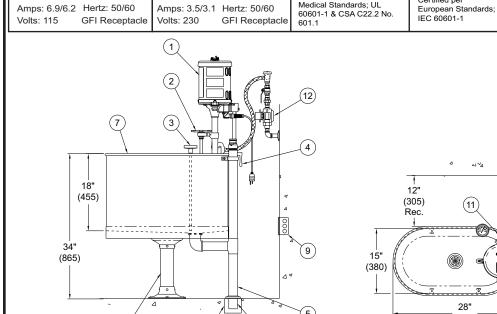

(10)

ELECTRIC REQUIREMENTS:

12" (305) (11)Rec. 28' (710)

NOTES:

- 1. 1/2" HP Jet Pump Motor
- 2. Variable Pressure Control
- 3. Thermometer
- 4. Turbine Raising & Lowering Device

STATIONARY E-22-SP SERIES WHIRLPOOL

Certified per US & Canadian

5. Combination Drain & Overflow Assembly

(8)

- 6. Mounting Pedestal
- 7. Stainless Steel Tank
- 8. 2" O.D. x 5" Tailpiece
- 9. 3 Pole GFI Receptacle by Others
- 10. Floor Drain by Others
- 11. Filler Spout
- 12. Mixing Valve, Shown For Reference Only

▲ WARNING: Cancer and Reproductive Harm - www.P65Warnings.ca.gov

All dimensions are nominal and subject to manufacturer's change without notice. Whitehall assumes no responsibility fo use of void or superseded data © Whitehall Manufacturing, Member of Morris Group International.

Please visit www.whitehallmfg.com for most current specifications

SELECTION SUMMARY & APPROVAL FOR MANUFACTURING		
Model Number & Options		Quantity
Company		Date
Contact	Title	
Approval for Manufacturing/Signature		

Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 782-7706 • (626) 968-6681 • Fax (626) 855-4862

Updated: 10/02/18