

Extremity Whirlpools

Division of Acorn Engineering Company



A MEMBER OF



Jse



E Series 45 Gallon Stationary Whirlpool

E-45-S



E-45-S

Whitehall E-45-S Stationary Whirlpool is recommended for treatment of the upper extremities and the lower leg below the knee.

Tank Capacity is 45 gallons.

Whirlpool Tank is 32"L x 15"W x 25"D and is seamless welded construction, fabricated from heavy gauge, type 304 stainless steel, polished to a satin finish. Tank rim is reinforced with a concealed stainless steel rod. Tank bottom is seamless, coved (rounded) design for effective cleaning to minimize bacteria buildup. Tank is also provided with an extended stainless steel skirt flush to the floor line to conceal internal drain piping.

A 2" Combination Drain and Overflow, and a Filler Spout are provided.

Turbine Raising and Lowering Device permits adjustment of desired height and direction of water agitation towards the area of the body to be treated. It is spring balanced and provided with a locking device.

Turbine Assembly s equipped with Underwriters Laboratory listed 1/2 HP jet pump motor protected with an ABS plastic shell, automatic thermal overload protector, and lifetime-sealed bearings. All parts contacting the water are stainless steel, chrome plated brass or bronze.

Thermometer is provided with a dual scale, adjustable calibration, and a stem inserts into a retainer to prevent rattle during whirlpool agitation.

GUIDE SPECIFICATIONS

Provide Whitehall Stationary Whirlpool (specify model number and options) with a 45-gallon capacity, 32"L x 15"W x 25"D. Whirlpool shall be fabricated from heavy gauge, type 304 stainless steel. Construction shall be seamless welded and exposed surfaces shall be polished to satin finish. Tank shall have a reinforced rim and shall be fitted with a stainless steel skirt flush to the floor line. Bottom of tank shall be a seamless, rounded coved design to minimize bacteria build up. Provide a 2" combination drain and overflow, and a filler spout. Provide a turbine assembly with raising and lowering device which functions as an agitator and permits adjustment of desired height, direction of water agitation, and can be locked into place. Turbine assembly shall be UL listed 1/2 HP jet pump motor with automatic thermal overload protector and lifetime-sealed bearings. Provide a thermometer with a stem retainer to prevent rattle when whirlpool is in operation.

Please visit www.whitehallmfg.com for most current specifications.





Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 782-7706 • (626) 968-6681 • Fax (626) 855-4862

Updated: 010/02/18



Extremity Whirlpools

MODEL NUMBER AND OPTION SELECTIONS

Series

☑ E Extremity Series

Capacity

₫ -45 Gallons

Type

Stationary Whirlpool

Turbine Electrical Requirements

- ☐ 115VAC / 50Hz ☐ 115VAC / 60HZ
- □ 230VAC / 50Hz □ 230VAC / 60HZ

Country of Use

For Mixing Valve:

Refer to supplemental MXT15 sheet for details.

Suffix Options (See Whirlpool Options Section for Details)

- ☐ -AHC1 Adjustable High Chair with Casters
- -AHC2 Adjustable High Chair with Rubber Tips

☐ -ARS Arm Rest Support

-EG Enviro Glaze Color Powder Coating

Specify Color:

(Special Finishes Consult Factory)

Ground Fault Interrupter

☐ -PCS Pump Strainer Cover☐ -RHG Rubber Hand Grips

☐ -SDO Separate Drain and Overflow Assembly

ELECTRIC REQUIREMENTS:


Amps: 6.9/6.2 Hertz: 50/60 Volts: 115 GFI Receptacle


ELECTRIC REQUIREMENTS: Amps: 3.5/3.1 Hertz: 50/60 Volts: 230 GFI Receptacle

Certified per US & Canadian Medical Standards; UL 60601-1 & CSA C22.2 No.

Certified per European Standards; **(E** IEC 60601-1



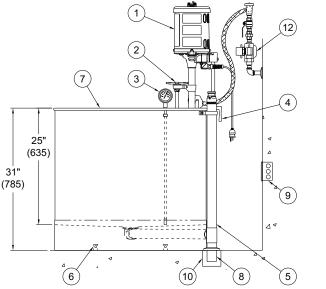

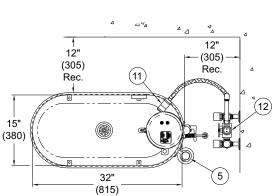

NOTES:

- 1. 1/2 HP Jet Pump Motor
- 2. Variable Pressure Control
- 3. Thermometer
- Turbine Raising & Lowering Device

STATIONARY E-45-S SERIES WHIRLPOOL

- Combination Drain & Overflow Assembly
- 6. "L" Brackets (4) Four Each
- 7. Stainless Steel Tank
- 8. 2" O.D. x 5" Tailpiece
- 9. 3 Pole GFI Receptacle by Others
- 10. Floor Drain by Others
- 11. Filler Spout
- 12. Mixing Valve, Shown For Reference Only

All dimensions are nominal and subject to manufacturer's change without notice. Whitehall assumes no responsibility fo use of void or superseded data © Whitehall Manufacturing, Member of Morris Group International.

Please visit **www.whitehallmfg.com** for most current specifications.

▲ WARNING: Cancer and Reproductive Harm - www.P65Warnings.ca.gov

SELECTION SUMMARY & APPROVAL FOR MANUFACTURING		
Model Number & Options		Quantity
Company		Date
Contact	Title	
Approval for Manufacturing/Signature		

Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 782-7706 • (626) 968-6681 • Fax (626) 855-4862