

### **Extremity Whirlpools**

Division of Acorn Engineering Company







Use



# E Series 22 Gallon Mobile Whirlpool

E-22-M, E-22-MU



E-22-MU

Whitehall E-22-M Mobile Whirlpools are recommended for treatment of the foot and ankle area, and when provided with a raised undercarriage, may also be used in treatment of the upper extremities.

Tank Capacity is 22 gallons.

**E-22-MU** is a Mobile Whirlpool with Undercarriage of heavy duty construction and fixed height which raises the tank for use in treatment of the upper extremities.

Whirlpool Tank is  $28^{\prime\prime}$ L x  $15^{\prime\prime}$ W x  $18^{\prime\prime}$ D and is seamless welded construction, fabricated from heavy gauge, type 304 stainless steel, polished to a satin finish. Tank rim is reinforced with a concealed stainless steel rod. Tank bottom is a seamless, coved (rounded) design for effective cleaning to minimize bacteria buildup.

**Heavy Duty Swivel Casters** support all mobile models. Rear pair of casters are lockable for stabilizing whirlpool during treatment.

**Turbine Raising and Lowering Device** permits adjustment of desired height and direction of water agitation towards the area of the body to be treated. It is spring balanced and provided with a locking device.

**Turbine Assembly** functions both as agitator and emptying device. It is equipped with Underwriters Laboratory listed 1/2 HP jet pump motor protected with an ABS plastic shell, automatic thermal overload protector, and lifetime-sealed bearings. All parts contacting the water are stainless steel, chrome plated brass or bronze.

**Drain option -SDP** has a bottom drain, leading to a separate, concealed water evacuation pump. A hose is attached to the side of the unit to facilitate the evacuation.

**Thermometer** is provided with a dual scale, adjustable calibration, and a stem inserts into a retainer to prevent rattle during whirlpool agitation.

#### **GUIDE SPECIFICATIONS**

Provide Whitehall Mobile Whirlpool (specify model number and options) with a 22-gallon capacity, 28"L x 15"W x 18"D. Whirlpool shall be fabricated from heavy gauge, type 304 stainless steel. Construction shall be seamless welded and exposed surfaces shall be polished to satin finish. Bottom of tank shall be a seamless, rounded coved design to minimize bacteria build up. Whirlpool shall have heavy duty swivel casters with locking device on rear pair of casters. Provide a turbine assembly with raising and lowering device which functions both as agitator and emptying device, permits adjustment of desired height, direction of water agitation, and can be locked into place. Turbine assembly shall be UL listed 1/2 HP jet pump motor with automatic thermal overload protector and lifetime-sealed bearings. Provide a thermometer with a stem retainer to prevent rattling when whirlpool is in operation.

Please visit **www.whitehallmfg.com** for most current specifications.





Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 782-7706 • (626) 968-6681 • Fax (626) 855-4862

Updated: 10/02/18



## **Extremity Whirlpools**

### **MODEL NUMBER AND OPTION SELECTIONS** Series

☑ E Extremity Series

**Capacity** 

✓ -22 Gallons

**Type** 

☐ -M Mobile Whirlpool

☐ -M-MU Mobile Whirlpool with Undercarriage

**Turbine Electrical Requirements** 

☐ 115VAC / 50Hz ☐ 115VAC / 60HZ ☐ 230VAC / 50Hz □ 230VAC / 60HZ

Country of Use\_

GFI Receptacle

#### For Mixing Valve:

Volts: 115

**Suffix Options** (See Whirlpool Options Section for Details)

☐ -ADH Auxiliary Drain Hose Assembly

☐ -AHC1 Adjustable High Chair with Casters ☐ -AHC2 Adjustable High Chair with Rubber Tips

☐ -ALC Adjustable Low Chair -ALS Adjustable Low Stool

-ARS Arm Rest Support Enviro Glaze Color Powder Coating --EG

Specify Color:\_\_\_\_\_\_(Special Finishes Consult Factory)

☐ -GFI Ground Fault Interrupter -HL Hydro Lift (Refer to -HL Page) -PCS Pump Strainer Cover

☐ -SDP Separate Drain Pump

Refer	to	supplement

**ELECTRIC REQUIREMENTS:** 


Amps: 6.9/6.2 Hertz: 50/60

tal MXT15 sheet for details.


**ELECTRIC REQUIREMENTS:** Certified per US & Canadian Medical Standards: UL Amps: 3.5/3.1 Hertz: 50/60 60601-1 & CSA C22.2 No GFI Receptacle

European Standards; IEC 60601-1



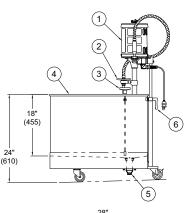


Volts: 230

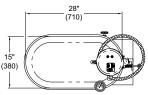

MODEL E-22-M

(2 (4) (3 (455) (6) (840)

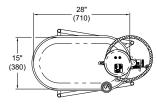

MODEL E-22-MU

### MOBILE E-22-M SERIES WHIRLPOOL

- 1. 1/2 hp Jet Pump Motor
- 2. Variable Pressure Control
- 3. Thermometer
- 4. Stainless Steel Tank
- 5. Auxiliary Drain Valve
- 6. Turbine Raising & Lowering Device
- 7. Optional -MU Under Carriage

MARNING: Cancer and Reproductive Harm - www.P65Warnings.ca.gov

All dimensions are nominal and subject to manufacturer's change without notice. Whitehall assumes no responsibility fo use of void or superseded data © Whitehall Manufacturing, Member of Morris Group International.

Please visit www.whitehallmfg.com for most current specifications

_	1			
SELECTION SUMMARY & APPROVAL FOR MANUFACTURING				
Model Number & Options	Quantity			
Company	Date			
Contact	Title			
Approval for Manufacturing/Signature				

Whitehall Mfg. • P.O. Box 3527 • City of Industry, CA 91744 • (800) 782-7706 • (626) 968-6681 • Fax (626) 855-4862